



ULTF-37i

ULT TECHNICAL SPECIFICATIONS



DIMENSIONS

Outer Dimensions HxWxD, mm	631x630x675
Depth incl. door handle	745
Inner Dimensions HxWxD, mm	433x439x271
Weight Gross/Net, kg	70 / 59
Material inner cabinet	Stainless steel
Material outer cabinet	Painted Steel
Packaging weight, kg	11
Packaging dimensions HxWxD, mm	850x800x850
Insulation thickness	80
Insulation type	Polyurethane with Cyclopentane
Mobility	Standard: Adjustable Feet / Option: Castors
Refrigerant, Type / gram	Nature R 2/96gr

FEATURES

Lock	☒
LED light	☒
Battery Backup for Controller, 24h	☒
Porthole	☒
Porthole size	20 mm
Dry contact	☒
Door reversibility	☐
Vacuum valve	☒
VIP (Vacuum Insulated Panel)	☒

STORAGE

Capacity (L)	37
--------------	----

ALARMS

High / Low temperature	☒
Open door	☒
Power failure	☒
Probe failure	☒

COOLING

Refrigerant	Nature R / 96
Noise	< 48
Voltage	100-240 V
Frequency	50/60 Hz
Temperature range	-20 / -86
Ambient condition	25 °C



ULTF-37i

ULT TECHNICAL SPECIFICATIONS




Voltage/Frequency	Voltage/Hz	100-240 V, 50/60Hz
Max Ambient	°C	32
Max Humidity	% rh	65
PERFORMANCE		
All data in RT20°C		
Temperature Range	°C	-20 to -86
Uniformity in performance - difference between top and bottom	°C	+/- 2,5
Pull down time	Minutes	84 min. to -75°C
Hold over time	Minutes	N/A
Noise	dB	< 48
Energy Consumption, kWh / 24h	kWh/24h	2,94 kWh/24h Set -80
Energy year	kWh/year	1186,6kWh/y Set -82°C
Instant Power Consumption	kW	PD 0,300/Stable 0,190
Heat Rejection*	W	150
U-Value	W/m ² K	N/A
COOLING COMPONENTS		
Refrigerant/Amount (gram)		Nature R 2/96gr
Number of compressors	pcs	1
Variable speed compressor	Yes/No	Yes
Internal air distribution (type of)		Static
Evaporator Fan	Yes/No/Variable	No
Condenser Fan	Yes/No/Variable	Yes
Number of probes	pcs	2
Defrost	Yes/No	No
CONTROLLER		
Controller	y/n/optional	yes
Controller type		i-Care touch
USB Connection	y/n	yes
Data connection	USB	yes
Controller abilities		logging of data & alarms
Controller languages		EN, DE, FR
Temperature curve in controller	y/n	yes

* Heat Rejection is defined as average power based on energy consumption, rounded up to nearest 50W.



ULTF-37i

ULT TECHNICAL SPECIFICATIONS

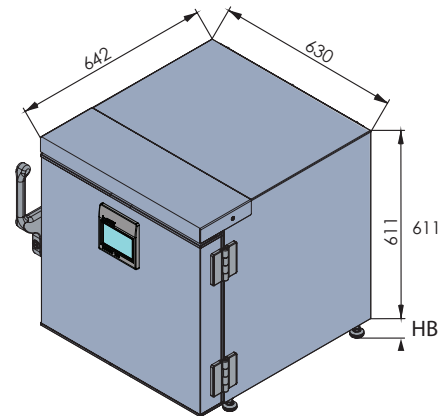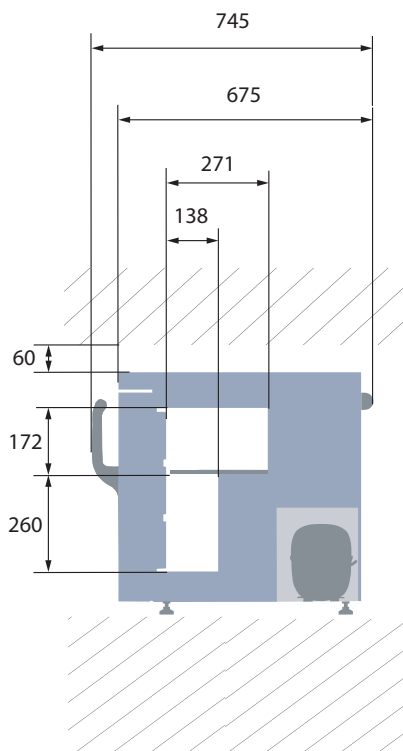
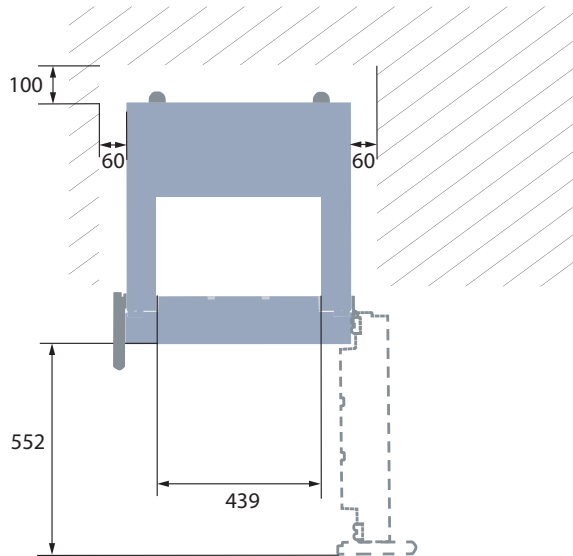


DIMENSIONS

Height, mm	631
Width, mm	630
Depth, mm	675
Depth incl. handle, mm	745

All measurements in mm.

HB: Height of base
(HB is adjustable when given value is xx-xx)

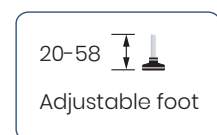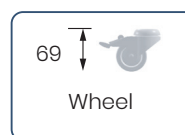




ULTF-37i

ULT TECHNICAL SPECIFICATIONS

SENSOR POSITION

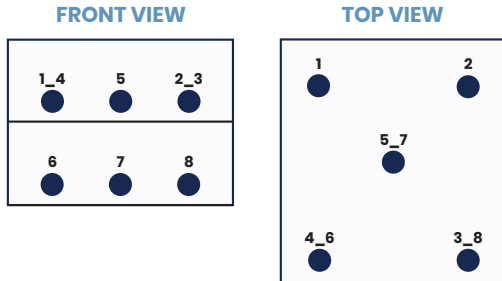

MODEL


ULTF-37i

Test type	8-point test
Test environment	Controlled conditions, empty cabinet
Ambient temperature	20°C
Humidity	60%
Set-point	-82°C
Sensor used	25gr tinned brass formed as a cylinder with a diameter of 15,2 mm
Installation	Appliance installed according to instruction manual conditions
Refrigerant	Nature R


SENSOR TEMPERATURE

Sensor position	P1	P2	P3	P4	P5	P6	P7	P8
Max	-79,7	-80,2	-80,1	-79,9	-80	-83,1	-83,1	-84,6
Avg.	-80	-80,5	-80,4	-80,2	-80,3	-83,5	-83,6	-85
Min.	-80,4	-81	-80,8	-80,5	-80,7	-84,1	-84,1	-85,4

WARM UP & PULL DOWN

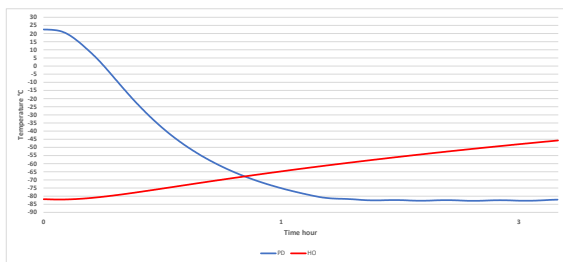

CYCLIC OPERATION

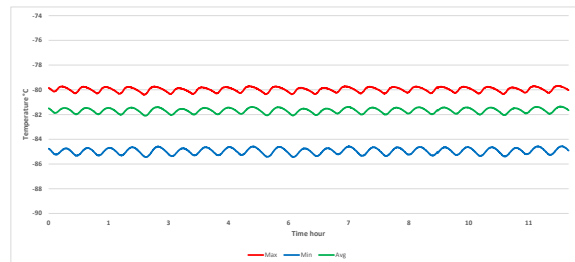

TYPICAL PERFORMANCE IN AMBIENT 20°C - EMPTY CABINET

Avg. cabinet temperature	-81,7°C
Peak variation from set-point	+/- 2,5°C
Stability in avg.	0,5°C
1 min. door open recovery to -75°C avg. temperature	20 min.
Cycle rate on/off	25,7 / 6,7 min.
Duty cycle	78,7%
Energy consumption - Normal mode	3,25 kWh/day
Energy consumption - Energy saving mode	2,92 kWh/day
Pull down time to -75°C avg. temperature	84 min.
Hold over time from -82°C to -60°C	105 min.
Heat rejection	306 W